
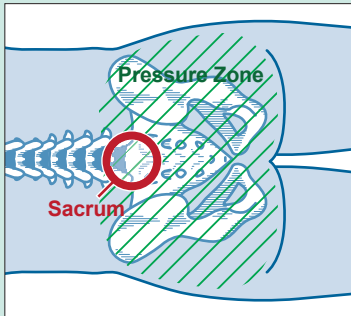
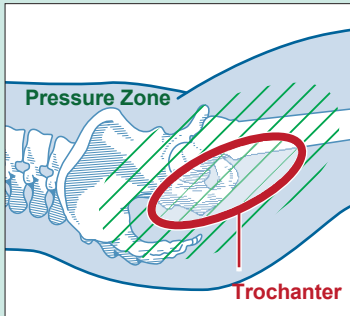
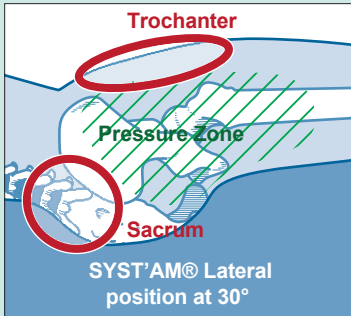


30° Positioning Wedge

Memory foam and a high resiliency stabilising foam nub






The 30° positioning wedge consists of memory foam and a high resiliency stabilising foam nub is widely recommended for those at risk of developing pressure ulcers due to laying in a supine position for an extended period of time. The wedge helps to protect the sacrum and trochanter by redistributing the user's body weight to lower risk areas of the body and also facilitates handling during repositioning.

Features	Benefits
<p>Viscoelastic foam with memory effect</p>	<p>Viscoelastic foam gives the positioning wedge better resistance against subsidence. The memory foam allows precise moulding which helps to provide effective pressure distribution and comfort for the user.</p> 
<p>The Form</p>	<p>The positioning wedge has lower limb support at thigh level and a hollowed sacrum zone to help redistribute pressure to lower risk areas of the user's body. The supports at thigh level also enhance user comfort by avoiding the rotation of the vertebral column, contributing to more effective use. It also prevents contact between the knees which can be at a high risk of pressure ulcers in the occurrence of prolonged contact.</p> <div style="display: flex; justify-content: space-around;"> <div data-bbox="331 1576 683 1890">  <p>Pressure on the sacrum</p> </div> <div data-bbox="740 1576 1091 1890">  <p>Pressure on the trochanter</p> </div> <div data-bbox="1145 1576 1497 1890">  <p>SYST'AM® Lateral position at 30° No pressure on the sacrum or trochanter</p> </div> </div>

UHC product code: USYC0WED1R (right) USYC0WED1L (left)
 SYST'AM® product code: P910LD1HW (right) P910LG1HW (left)

Stabilisation Insert

The 30° positioning wedge has a central foam anatomical nub which is moulded to a soft viscoelastic surface. The nub is a key feature for the effective use of the wedge as it enables the 30° position to be maintained over time by acting as a support structure. The ergonomic shape of the insert acknowledges the curves of the user's body which helps to reduce pressure.

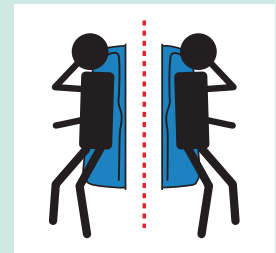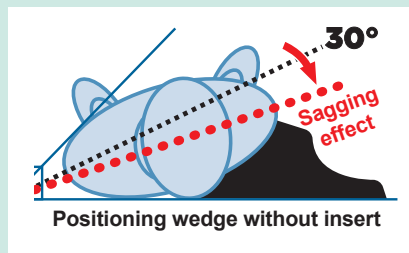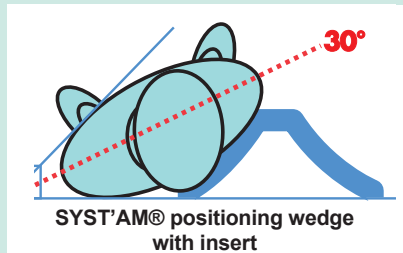

The Cover

The breathable, waterproof and vapour permeable cover promotes comfort whilst aiding infection control. The wedge is covered with a soft, bi-stretch polyurethane cover and features a non-slip base and carry handle to facilitate transfers from one room to another. The cover is also marked with a marker label to ensure perfect positioning.



Technical Specification

UHC Product code	SYST'AM® Product code	Cushion Dimensions	Note
USYC0WED1R (right wedge)	P910LD1HW	55 x 25 x 7.5cm / 21½ x 10 x 3"	For a user turning to the left
USYC0WED1L (left wedge)	P910LG1HW	55 x 25 x 7.5cm / 21½ x 10 x 3"	For a user turning to the right

Warranty: 2 years

Cleaning of the cover:



Cleaning of foam:

Do not immerse in water



UHC product code: USYC0WED1R (right) USYC0WED1L (left)
SYST'AM® product code: P910LD1HW (right) P910LG1HW (left)