

Diversi-foam Hybrid Cushion

Dual combination self adjusting cushion


Coccyx cutout

Coccyx air/ foam cell

Ischial air/ foam cell

Rear thigh air/ foam cell


Front thigh air/ foam cell

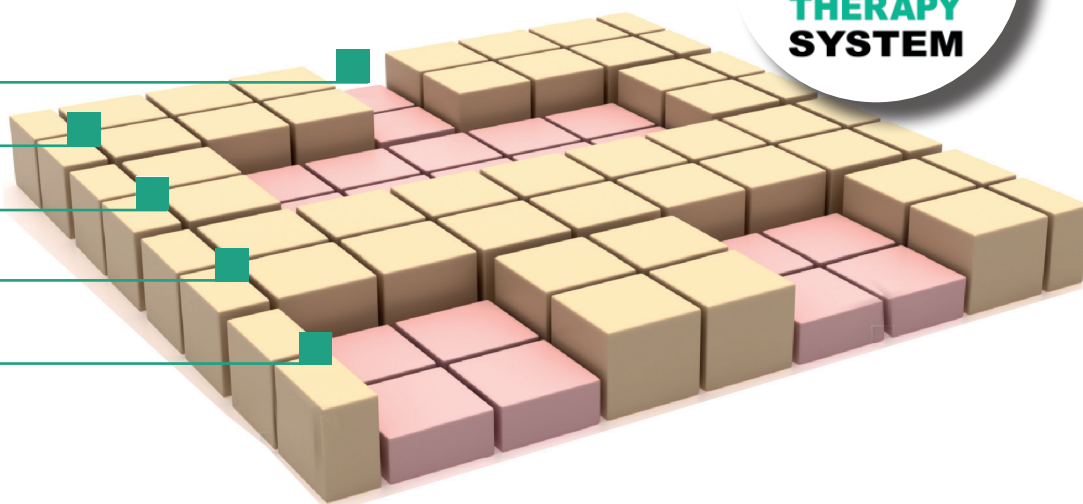

The **Diversi-foam hybrid cushion** incorporates air and foam filled cells with visco elastic inserts to redistribute pressure away from high risk areas without requiring initial set up or regular intervention.

The cushion adjusts automatically to the users weight distribution and can be used for patients up to a **High Risk** of pressure ulcer development, within multiple care settings.

Features & Benefits

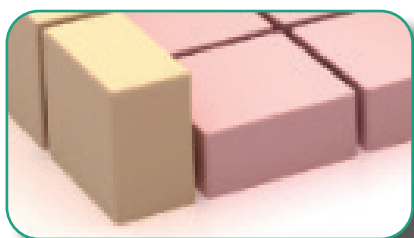

Foam inserts

Visco elastic foam allows for immersion, comfort and postural support whilst re-distributing pressure away from high risk areas



Coccyx, Ischial, front and rear thigh air cells

Castellated visco elastic foam within aircells offering continuous seated support and pressure redistribution



Ultimate Immersion therapy

A combination of aircells and an intrusive valve react to provide natural user immersion. The negative pressure created also benefits the user with additional postural support

For more information on this product please contact us on
tel: 0333 321 8996 or email: sales@ultimatehealthcare.co.uk

www.ultimatehealthcare.co.uk





The **Cushion Cover** is waterproof and vapour permeable with an enclosed zip to maintain high infection control standards.

The Diversi-foam Cushion has the option for 24 hour care by using in combination with the Diversi-foam II Hybrid Mattress system

Specification

Available in multiple sizes and depths. Please see size chart below

Product code:	UPRH1717-2	Dimensions:	43cm x 43cm x 5cm	(17" x 17" x 2")
	UPRH1717-3		43cm x 43cm x 7.5cm	(17" x 17" x 3")
	UPRH2020-2		50 cm x 50 cm x 5 cm	(20" x 20" x 2")
	UPRH2020-3		50 cm x 50 cm x 7.5 cm	(20" x 20" x 3")
Risk category:	Up to a high risk	Warranty:	2 years	
Max user weight:	152kg 24 stone 160kg 25 stone	Product weight:	0.9kg - 1.4kg	

Cleaning: Can withstand up to 10,000 ppm chlorine (to be used in conjunction with local infection control guidelines) (can be laundered up to 84°C)



For more information on this product please contact us on
tel: 0333 321 8996 or email: sales@ultimatehealthcare.co.uk

www.ultimatehealthcare.co.uk