

## In-line Grab Rail




The In-line Grab Rail provides a stable handle for patients to hold while entering or exiting the bed.

The In-line Grab Rail can be placed anywhere along the bed frame however there must be gaps of less than 60mm or more than 318mm from the head and foot boards to reduce the potential risk of entrapment.

The In-line Grab Rail is suitable for patients with some mobility who require additional support whilst standing up from the bed or sitting back down. It can be secured using 2 easy grip hand knobs.

**If risk assessed can be used in conjunction with  
Full length side rails & 3/4 length side rails**

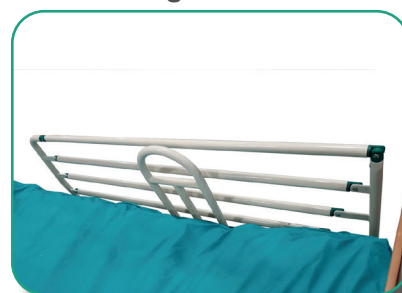

The in-line grab rail is only suitable for use on Cura II, Cura II Extra High, Cura Supreme and Cura Community Beds

### Features

- Fixes securely to bed frame
- Helps maintain independence
- Multiple positioning
- Facilitates entry and exiting of bed

### Specifications

<b>Product code:</b>	UBEB7004-3	<b>Overall height (from platform):</b>	450mm
<b>Overall height:</b>	628mm	<b>Overall width:</b>	225mm

For more information on this product please contact us on  
tel: 0333 321 8996 or email: [sales@ultimatehealthcare.co.uk](mailto:sales@ultimatehealthcare.co.uk)

[www.ultimatehealthcare.co.uk](http://www.ultimatehealthcare.co.uk)

