


Heritage II Digital Turn

Click to play video

Dynamic Mattress Replacement Incorporating Continuous Turning and Turn Assist

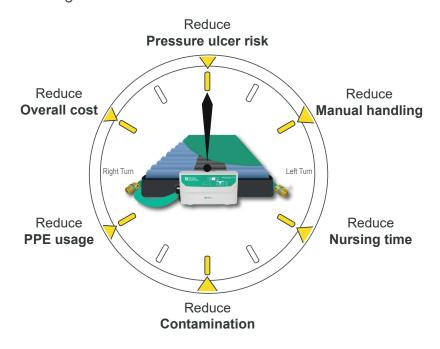

The Heritage II Digital Turn Dynamic Mattress Replacement system provides very high optimal therapy for the prevention and treatment of pressure ulcers in any care environment.

Manual turning of patients can be intrusive, especially during the night and often results in disturbed sleep.

The Heritage II Digital Turn mattress can reduce the care package required.

The mattress has the option of right and left or bilateral turn to assist in a manual or provide automatic turning at variable times between 10 - 30 minutes.

Traditional Mattress Systems



Require double handed care for the turning of patients which can be intrusive, disturb sleep and result in carer and patient injury.

Single Handed Care Turn Assist



The amount of effort needed to manually turn a patient is significantly reduced, taking away the need for a second carer.

Features & Benefits

Continuous



Reducing the need for care, where a patient needs to be turned on a regular basis.


Activation of the Turning Valve



Allows for both assisted turning when required or automatic unassisted continuous turning reducing care support.



Low Air Loss - Climate control

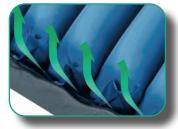

Reducing skin heat and moisture by utilising the alternating air after each cycle releasing up to 100 litres of air underneath the cover to ventilate and cool the skin.

Power Unit

- Digital control panel
- Comfort control
- Auto-Detection pressure setting Detecting the patients weight distribution to set internal pressures appropriately.
- Dual 9 LPM compressors
- Alternating 2-cell cycle
- Variable cycle times (10/15/20/25/30 mins)
- Automatic control panel lockout
- Timed static mode (20 minutes)

Automatically returning to alternating mode.

Auto-Firm (20 minutes)

Automatically retuning to the Last Mode of Operation

- Transport mode Up to 24 hours
- · Audible and visual alarms
- · Annual service indicator
- · Adjustable hanging brackets
- · Constant low pressure mode -

Allows the mattress to conform closely a patient's body and redistribute body weight.

TruComfort™



A 1 in 2 alternating system, where the deflated cell fully re-inflates before the next cell deflates. Eliminating the movement some patients experience with the traditional alternating cells.

Mattress

- · Left and right turning
- · Cell-on-cell construction -

Constructed of deep 8" cells, providing dual therapy modes i.e. alternating and constant low pressure

- Static base and head cells Providing support and stability
- Lower leg zone
- CPR Valve
- Strong cell retention knitted bands
- Maximum patient weight 254kg / 40 stone

Cover

- · Soft, waterproof vapour permeable PU cover
- · 360 degree zipper
- · Welded seams and enclosed zip
- · Integral covered air pipe
- Cable management system

Specification

Product code:

Risk category:
Up to very high risk

Max user weight:
254kg | 40 stone

Mattress dimensions:
2000 x 890 x 200mm

Power unit dimensions:
360 x 130 x 230mm

7.4kg mattress | 4.6kg pump

Cell material:Nylon PUCover material:Polyester PU

Cycle time: 10 /15 / 20 / 25 / 30 mins

(2-cell cycle)

Pressure range: 15-60mmHG Supply voltage: 230V, 50 Hz

Fuse rating: T1A

Current rating: 0.12 A

Classification: Class II, Type BF

Compliance: EN 60601-1, EN 60601-1-2

Warranty: 2 years

Pressure Mapped Patient



During right turn mode

Age: 44

Gender: Male Weight: 83kg Height: 170cm

For more information on this product please contact us on tel: 0333 321 8996 or email: sales@ultimatehealthcare.co.uk

